

# Focus Acrylic™

*By  
eyePx*



eyePx - Model 601 Hydrophobic Acrylic  
Aspheric Foldable Intraocular Lens  
UV Blocker

# Focus Acrylic™


## MODEL 601

*An acrylic IOL designed to give you all the advantages of hydrophobic acrylic IOL's with none of the compromises.  
We invite you to evaluate Model 601*


### Focus Acrylic One Piece Aspheric Hydrophobic Lens

Model 601 is designed as a Single Piece Lens of highly biocompatible hydrophobic acrylic material with advanced optics, and UV blocker. The haptic also manufactured from the same acrylic material with 0 degrees angulation and a total diameter of 13.0mm with a full 6.0mm square edge biconvex aspheric optic. The color of the lens is transparent and it is injectable through a 2.2mm incision using the titanium EL-22M Epsilon Injector in conjunction with the OD522 Ophtec cartridge. For the disposable option, use the Ophtec OD665 cartridge and injector or Comport P RET Inc. The lens is designed to be positioned in the posterior chamber. Available diopters from 0.0 to 34.0

#### Light Transmission Curves

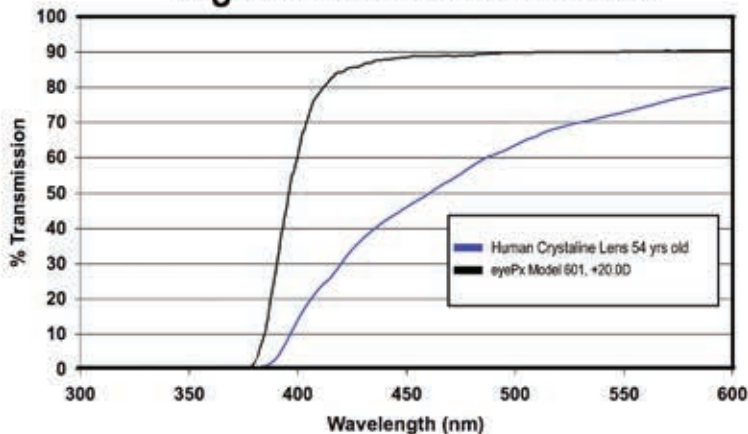

#### Proven Design Properties

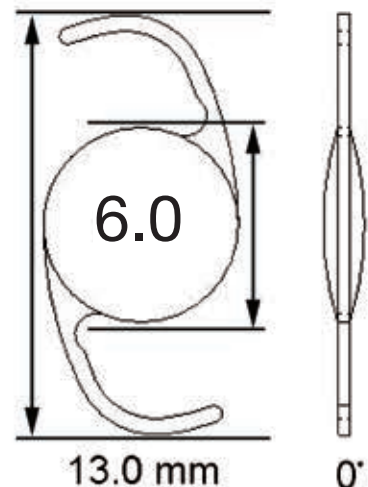

Optic: Full symmetric biconvex aspheric  
Material: Hydrophobic acrylic  
Range: 0.0 to 34.0 Diopters  
Haptic: Acrylic — Modified L  
A Constant: 118.3 (estimated)  
IOL Master: 119.0 (calculated)  
Anterior Chamber Depth: 5mm



Model 601 may be implanted with an Epsilon EL-22M injector

eyePx

Manufactured by eyePx, LLC. 92340 Bourg La Reine, France  
email: [eyePx@cox.net](mailto:eyePx@cox.net)

CE  
PIN 100489-001  
REV. C